



## USER GUIDE



**ATTENTION:** Please read this guide before using your hearing aids. Following these steps will help you achieve the best results.



# HEARING AID

MODEL: JH-490  
ART NO.: JH-490

**MD** OTC

**WARNING:** This product is not intended for individuals under 18 years old. If you are younger than 18, please consult a doctor, ideally an ear, nose, and throat specialist (ENT), as your condition requires professional medical attention. Over-the-counter hearing aids are designed exclusively for adults aged 18 and older.

This OTC hearing aid is intended for users aged 18 and older. If you are under 18 and experiencing hearing loss, it is important to consult a doctor, preferably an ear, nose, and throat specialist (ENT). Younger individuals may require medical testing and specialized management. Since hearing loss can impact speech and learning, professional fitting and ongoing care are essential.

## **WARNING: When to Consult a Doctor**

If you experience any of the following issues, seek medical advice from a doctor, ideally an ENT specialist:

- Presence of blood, pus, or fluid from your ear within the past six months.
- Persistent ear pain or discomfort.

- Excessive earwax or foreign objects: If you suspect a lot of earwax or something in your ear, consult a doctor.
- Dizziness or vertigo: Seek help if you feel dizzy or experience a spinning sensation (known as vertigo).
- Sudden hearing changes: Get medical attention if your hearing has changed suddenly within the past six months.
- Fluctuating hearing loss: If your hearing worsens and then improves, it may require professional evaluation.
- Uneven hearing: Hearing significantly worse in one ear could be a sign of a condition needing medical attention.
- Ringing or buzzing in one ear: If tinnitus occurs only in one ear, consult a specialist.
- Ear deformities or injuries: For ears with birth defects, unusual shapes, or injuries, professional advice is recommended.

**WARNING: Pain or discomfort while inserting the hearing aid**

- Remove the device immediately if it causes pain or discomfort when inserting. Follow the provided instructions before trying again. If pain or discomfort persists, contact the manufacturer or a hearing healthcare professional. You can report this as an adverse event to the FDA using the instructions included.

**Caution: This device is not a substitute for hearing protection.**

If you are in a loud environment, such as a concert or construction site, remove the device and use proper ear protection like earplugs. Prolonged exposure to loud sounds while wearing this device can be harmful.

**Caution: Avoid painful or uncomfortable sound output.**

Reduce the volume or remove the device if the sound output feels uncomfortable or painful. Persistent issues with volume adjustments may indicate the device needs further calibration.

**Caution: Medical assistance may be required if a piece gets stuck.**

If any part of the device, such as the ear-tip, becomes lodged in your ear and cannot be easily removed, seek medical help immediately. Do not use tweezers or cotton swabs, as they could push the part deeper into your ear, potentially causing serious injury.

**Note: Consult a professional if concerns persist.**

If you experience ongoing difficulties or concerns with this device, contact a hearing healthcare professional for evaluation and assistance.

If you experience hearing loss in both ears, using hearing aids in both ears can significantly enhance your listening experience, especially in challenging environments like noisy settings that can cause listening fatigue.

**Note: Reporting Injuries, Malfunctions, or Adverse Events to the FDA**

If you encounter a problem with your hearing aid, report it to the FDA as soon as possible. These issues, referred to as "adverse events," may include:

- Skin irritation in or around the ear.
- Injuries caused by the device, such as cuts, scratches, or burns from overheating batteries.
- Device parts becoming lodged in the ear.
- Sudden worsening of hearing loss after using the device.
- For instructions on how to report, visit FDA MedWatch, call 1-800-FDA-1088, or download a form to submit by mail.

# TABLE OF CONTENTS

## (1) **RETRAINING YOUR BRAIN**

(5) Week 1

(10) Week 2

(12) Week 3

(14) Week 4

## (17) **QUICK START GUIDE**

(20) In the Box

(22) Meet Your Hearing Aids

- (24) Identifying Left and Right  
Hearing Aids
- (26) Charging Your Case and Hearing  
Aids
- (30) Understanding the Case Lights
- (39) HOW TO WEAR  
YOUR HEARING AIDS
- (44) Turning On and Adjusting Your  
Hearing Aids
- (48) Adjusting the Volume  
Changing Modes



- (51) Selecting the Right Ear Dome
- (59) Removing Your Hearing Aids
- (61) Protecting Your Hearing Aids
- (63) **TROUBLESHOOTING**
- (87) Care and Cleaning

# ADJUSTING TO YOUR HEARING AIDS:

---

Your 30-day plan for the best  
hearing experience possible.

# HOW TO GET STARTED

The Oricle 3.0 is an advanced hearing device designed to enhance your hearing. However, your hearing may need time to adjust.



To ensure you get the best results, we've created a detailed 30-day adjustment plan tailored for all our customers.

# PLEASE FOLLOW THESE STEPS DURING THE FIRST 30 DAYS



# BEFORE GETTING STARTED

**Prepare Your Ears:** Ensure your ear canal is free from wax buildup before using your new device. Consider a professional ear cleaning if necessary.

**Charge Your Hearing Aids:** Make sure your hearing aids are fully charged before their first use. (Refer to page 26 for detailed charging instructions).



## Day 1

- Start Slow: Wear your hearing aids for 1-2 hours daily.
- Adjust volume gradually for comfort.
- Start in quiet settings at home and avoid noisy environments.
- Practice one-on-one conversations in quiet environments.



## PLEASE NOTE:

**On the first day, you may notice sounds you haven't heard in a while**, like the hum of the refrigerator, the ticking of a clock, or the air conditioner blowing.

This might feel a bit distracting initially, but it's completely normal. These sounds are "new" to you, so your attention may naturally focus on them more than usual.

Over time, your brain will naturally filter out these sounds.

- **Feedback/Chirping:** It is normal to experience feedback or chirping when the hearing aid is inserted into your ear while powered on.
- Feedback may also occur if your hand is placed near your ear while the hearing aid is active.





## Days 2 - 7 (Week 1):

- Gradually increase the number of hours you wear your hearing aids each day.
- Continue using them only in quiet environments to help your hearing adapt to these "new" sounds that you may not have heard in a long time.

- Encourage your friends and family to support you by speaking clearly and at a steady pace. It's also helpful if they face you while talking, allowing you to align their mouth movements with the sounds and words you hear.
- **Adapting to hearing aids is a unique journey for everyone. Stay patient and give yourself time to adjust as you become more comfortable with them.**



## Days 8 - 14 (Week 2):

- As you continue wearing your hearing aids for several hours each day, begin exploring a variety of sounds, such as watching TV or enjoying music.
- Challenge yourself by listening to an audiobook or podcast to further enhance your auditory adjustment.

- Try watching movies or TV without relying on subtitles to strengthen your listening skills.
- This week, aim to wear your Oricle 3.0 hearing aids as often as possible to help your hearing adapt and sharpen its focus on key sounds at these enhanced levels.



## **Days 15 - 21 (Week 3):**

- In your third week, begin using your Oricle 3.0 hearing aids outdoors. By now, you are becoming accustomed to their benefits.
- Take walks and chat with people, but avoid noisy places like restaurants as they may take longer to adjust to.

## Quick Challenge:

Have a one-minute conversation with your hearing aids inserted, then remove them and continue talking. Compare the difference!



Has wearing your hearing aids started to feel natural and effortless?



## Days 22 - 28 (Week 4):

- After three weeks of consistent use, start venturing into more challenging environments like restaurants and supermarkets.
- Remember, noisy settings can be difficult not just for hearing aid users but even for those with normal hearing. Take your time adjusting.

Be patient as you explore what you can hear and understand in noisy places. It's normal to encounter some challenging or confusing sounds.





## IMPORTANT REMINDER:

We understand your eagerness to get back to normal hearing quickly. **However, jumping straight to Week 4 without following the recommended steps can lead to two outcomes:**

1. You may feel like your hearing aid isn't functioning properly.
2. Your ears will require significantly more time to adapt to the device.

# QUICK START GUIDE

---

## **KEY INFO**

**Inside the next pages, you'll find a Quick Start Guide to help you get started.**

**Access it here:**

**[support@oriclehearing.com](mailto:support@oriclehearing.com)**

**Or Scan the QR code**




# IN THE BOX



**Hearing Aids**



**Charging Case**

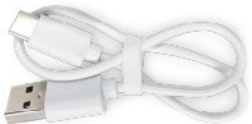

**Charging Cable**



**Ear-tips**



**Hearing Aid Brush**

# MEET YOUR HEARING AIDS

1. Microphone
2. Faceplate
3. Removal Handle
4. Charging Contact
5. Multi-Function Button
6. Ear-tip
7. Ear-dome



Multi-Function  
Button





# HOW TO TELL LEFT FROM RIGHT



**LEFT**

# HOW TO TELL LEFT FROM RIGHT



**RIGHT**

## **STEP 1: CHARGE YOUR CASE AND HEARING AIDS**

Before first use, ensure your hearing aids are fully charged.

### **IMPORTANT:**

Your hearing aids do not come with an adapter. You can use standard adapter from home (such as your phone charger) or connect the charging cable directly to a power bank or computer.

**1.) Connect the USB end of the cable to a wall adapter.**



**2.) Make sure the case and hearing aids are connected to charge.**



### 3.) Place the hearing aids in the case.



# WHAT DO THE 3 CASE LIGHTS MEAN?

- **Left Light:** Indicates the left hearing aid's battery level.



- **Right Light:** Indicates the right hearing aid's battery level.





- **Green Light:** Fully charged.



- **Orange Light:** Charging.



# EXAMPLES

LEFT

RIGHT

|                                                                                   |                  |                  |
|-----------------------------------------------------------------------------------|------------------|------------------|
|  | Fully<br>Charged | Fully<br>Charged |
|  | Charging         | Fully<br>Charged |
|  | Charging         | Fully<br>Charged |
|  | Charging         | Charging         |

# EXAMPLES

|                                                                                                                                                                     | LEFT             | RIGHT    |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------|----------|
|   | Fully<br>Charged | Charging |
|   | Fully<br>Charged | Charging |

- Ensure the hearing aids are correctly positioned in the charging case.



- They will magnetically click into place when properly aligned.

- A fully charged case provides 2-3 charging cycles for a pair of hearing aids.
- If unplugged, case lights will turn off 5 minutes after the hearing aids are fully charged.

- **Low Battery Indicator:** If the middle light is flashing orange the battery is <30% and needs to be charged.



- Hearing Aid Usage Time by Charging Duration:

**Charging Time**

**Expected Hearing  
Aid Usage Time**

|         |             |
|---------|-------------|
| 2 Hours | 10-14 Hours |
|---------|-------------|

## **STEP 2: HOW TO WEAR YOUR HEARING AIDS**

- **Insert the hearing aid**
- Make sure the right aid (red) goes in the right ear and the left aid (blue) in the left ear (*refer to Page 24*).
- Insert hearing aid with the removal cord pointing downward.





**Left**

**Right**

## **KEY INFO**


- Gently insert the hearing aid into the ear canal.  
To prevent the hearing aid from sitting too close to the eardrum, do not push the faceplate deeper than the opening of the ear canal.

- Use your finger or thumb to gently position the hearing aid in the ear if needed. The hearing aid may need to be slightly rotated once in the ear to obtain correct placement.

- If you experience feedback or chirping during insertion, turn the hearing aid off by pressing the button for 2 seconds.
- Once inserted, turn the hearing aid back on.

## **STEP 3: TURNING ON AND ADJUSTING YOUR HEARING AIDS**

Learn how the multi-function button works:

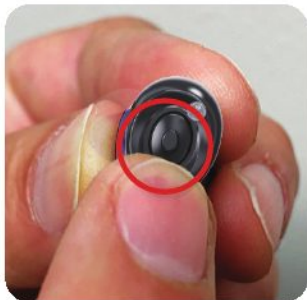

### **TURN ON AND OFF**

Press and Hold

### **VOLUME CHANGE**

Quickly Touch

### **PROGRAM CHANGE**

Quickly Press Twice 2X

## 1.) Powering Your Hearing Aid ON and OFF

- To turn the hearing aid **ON**: Press and hold the button for **2 seconds** until it powers on.
- To turn the hearing aid **OFF**: Press and hold the button for **2 seconds** until it powers off.
- Hearing aids automatically turn off when put them into the charging case.



## 2.) Adjusting the Volume

- **Tap the button** briefly to increase the volume, cycling through **levels 1 to 5.**







## PLEASE NOTE:

- **1 Du:** Indicates the volume has increased by one level.
- **1 Beep:** Signals the maximum volume (level 5).
- After reaching max volume, the next tap will reset the volume to level 1.

### 3.) Changing the Mode

- The Oricle 3.0 offers 2 modes.  
To switch modes, press the button quickly twice.

- The number of beeps will indicate the current mode.



## Modes and Their Beeps:

**1 Beep – Normal Mode:** Ideal for quiet spaces like home or office.

**2 Beeps – Restaurant Mode:** Best for noisy settings like restaurants or meetings.

 **KEY INFO**

The hearing aids are not synced, allowing independent adjustments for each ear.

**Use the button on each device separately to customize the settings to your preference.**

## **STEP 4: CHOOSE THE RIGHT EAR DOMES**

Now that you're familiar with using your hearing aids, it's time to select the ideal dome.

**Keep in mind that the dome can impact:**

- Comfort
- How your voice sounds
- The overall sound amplification
- Feedback or ringing

Finding the right dome may require some trial and error to match your ears and listening preferences.



## 1.) Attach the Dome:

Firmly attach a new dome to the tip of the hearing aid as illustrated. Gently tug it afterward to ensure it is securely in place.



**2.) With the dome attached, place the powered-on hearing aid into your ear:**



### **3.) Ensure a Secure Fit:**

Open and close your mouth a few times to confirm the hearing aid stays in place.





 **KEY INFO**

If the dome is not properly secured before use, it may detach and become stuck in your ear canal when removing the hearing aid. If this occurs, we advise seeking professional assistance for safe removal.

**Note: Ear domes need to be changed monthly and cleaned weekly. For best performance, clean them daily.**

# **DOME TROUBLESHOOTING TIPS**

## **Common Issues    Suggestions**

|                                                                 |                             |
|-----------------------------------------------------------------|-----------------------------|
| <b>Hearing aid<br/>keeps falling out<br/>of the ear.</b>        | Try a bigger<br>dome        |
| <b>Hearing aid<br/>won't fit; the<br/>dome is too<br/>large</b> | Switch to a<br>smaller dome |

# **DOME TROUBLESHOOTING TIPS**

## **Common Issues    Suggestions**

|                                                   |                                                                    |
|---------------------------------------------------|--------------------------------------------------------------------|
| <b>Feedback or chirping from the hearing aid.</b> | Ensure a snug fit or try a different dome size or style.           |
| <b>Your own voice sounds uncomfortable</b>        | Change the dome size or type for better comfort and clearer sound. |

## **STEP 5: REMOVING YOUR HEARING AIDS**

Hold the removal handle and gently pull the hearing aid out.



## **KEY INFO**

After removing your hearing aids, please put them in the charging case when not in use, and they will automatically turn off.



## How to Protect Your Hearing Aids

- **Keep Dry:** Remove hearing aids before showering, swimming, or any activity involving water or moisture.
- **Avoid Sprays:** Take them off before applying perfume or hairspray.

- **Store Safely:** Keep them out of reach of pets and children, and always store them in the charging case when not in use.
- **Temperature Matters:** Avoid exposing them to extreme heat, cold, or humid conditions.
- **If Wet:** Let them air dry. Do not use a blow dryer, microwave, or heat source.

# TROUBLESHOOTING

---



# TABLE OF CONTENTS

- (65) Feedback / Chirping
- (71) Barrel-like sound
- (75) Low volume
- (79) No sound
- (81) Charging issues
- (83) Chirping during insertion
- (84) Chirping when using the phone
- (85) Background noise
- (87) Maintenance and cleaning
- (90) Replacing parts

# TROUBLESHOOTING: FEEDBACK/CHIRPING

## **Why It Happens:**

Feedback or chirping occurs when the amplified sound from the hearing aid is picked up by the microphone and amplified again. This creates a sound loop, leading to the high-pitched noise.

## COMMON CAUSES:

- 1. Poor Ear Tip Fit:** An improperly seated hearing aid allows sound to escape and re-enter the microphone.
- 2. Volume Too High:** Excessive volume causes sound to leak and be re-amplified.
- 3. Nearby Objects:** Hands, phones, or hats near the hearing aid can reflect sound back into the microphone.

## How to Fix Feedback/Chirping:

### **FIX #1: Secure the hearing aid properly in your ear.**

- Lift the top of your ear gently with one hand while inserting the hearing aid with the other. This helps align the ear canal for a snug fit.
- Adjust or rotate the hearing aid slightly in your ear to achieve a tighter seal and reduce sound leakage.

## **FIX #2 - Clean the Hearing Aid and Dome.**

- Detach the dome and check the hearing aid tip and dome for any wax or debris.
- Use a hearing aid brush or a dry cloth to gently remove any buildup.



### **FIX #3 - Clean the Microphone.**

- Use a hearing aid brush to gently clear wax or debris from the microphone on the faceplate.



**FIX #4 - Lower the Volume.**

**FIX #5 - Switch Dome Size or Type.**



# TROUBLESHOOTING: BARREL-LIKE SOUND

## **Why It Happens:**

The "barrel" effect occurs when the hearing aid obstructs the ear canal, causing your own voice or specific sounds to become trapped and overly amplified inside your ear.



## COMMON CAUSES:

- 1. Blocked Ear Canal:** A completely sealed ear canal traps natural vibrations from your voice or other sounds, reflecting them back to the eardrum.
- 2. Poor Fit:** An ill-fitting dome that seals the ear too tightly can intensify the barrel effect.

## How to Fix Barrel-like Sound

### **FIX #1: Reposition the Hearing Aid.**

- Adjust the fit by slightly moving the hearing aids outward if they feel too tight.

### **FIX #2: Change the Program Setting**

- Switch to a different program mode (refer to Page 48 for instructions).

**FIX #3: Use a Smaller Dome.**



**FIX #4: Remove the dome to reduce the barrel effect:**



# TROUBLESHOOTING: LOW VOLUME

## **Why It Happens:**

The volume might be low because it hasn't been adjusted to the maximum level, or wax and debris may be obstructing the sound port.

## How to Fix Low Volume:

### **Fix #1: Increase the Volume:**

- Ensure the volume is set to the maximum level (refer to Page 48 for instructions).

### **Fix #2: Change the Program:**

- Switch to a different program setting for better amplification.

### **Fix #3: Clean the Hearing Aid:**

- Remove the dome and check for wax or debris blocking the sound port or microphone. Use the included brush or a damp cloth to clean the dome, tip, and microphone.

## **Fix #4: Clear Ear Wax:**

- If the hearing aid whistles when cupped in your hand, it's functioning, but wax may be blocking your ear canal. Use an ear rinse kit or seek professional cleaning to remove wax buildup, which can worsen hearing loss and reduce the aid's effectiveness.

# TROUBLESHOOTING: NO SOUND

## How to Fix No Sound:

### **Fix #1:** Check Charging Status:

- Ensure the hearing aids are fully charged. Look for green lights on the charging case (refer to Page 26 for guidance).



## **Fix #2: Power On the Device:**

Hearing aids automatically turn on 12 seconds after removal from the charging case. (see Page 44-45 for details)

## **Fix #3: Clean the Device:**

- Remove the dome and inspect the sound port and microphone for wax or debris. Use the included brush or a damp cloth to clean the tip, dome, and microphone thoroughly.

## **TROUBLESHOOTING: CHARGING ISSUES**

1. Hearing aids show as charged (green lights) but only last 15 minutes.
2. Charging case shows orange lights, indicating the hearing aids didn't charge after being in the case for a long time.

## How to Fix Charging Issues:

**Fix #1:** Turn off the hearing aids before placing them in the charging case. Hearing aids left on will not charge properly (see Page 26 for detailed charging instructions).

**Fix #2:** Clean the charging contacts with a damp cloth and try again. Avoid using water to clean the copper contacts.

## TROUBLESHOOTING: CHIRPING DURING INSERTION

**Fix #1:** Turn the volume to a low setting before inserting the hearing aid. This prevents feedback during insertion. Once the hearing aid is properly in place, turn it back on.

**Fix #2:** Start on low volume.

**Fix #3:** Try different ear tip sizes until you find a snug fit without any feedback or whistling.

## TROUBLESHOOTING: CHIRPING WHEN USING THE PHONE

**Fix #1:** Hold the phone as usual. If feedback occurs, adjust by moving the phone's speaker slightly away from your ear or the hearing aid's microphone until you find the best position.

## **TROUBLESHOOTING: BACKGROUND NOISE**

### **Why It Happens:**

Hearing aids amplify all sounds, including faint environmental noises like a fan or air conditioner that you may not have noticed before. For those with hearing loss, these sounds might feel unfamiliar at first.

This adjustment period is normal as your brain learns to recognize amplified sounds.

**Fix #1:** Reduce the volume to make the noise less noticeable.

# **MAINTENANCE AND CLEANING:**

---



# MAINTENANCE AND CLEANING

## Daily Care:

Gently clean the dome, wax guards, buttons, and microphone using the provided brush to remove earwax, dust, or other debris.



Wipe the hearing aid body with a dry, soft cloth.  
Avoid using water, alcohol, or any liquid  
cleaners, as they can cause damage.



## REPLACING PARTS

**Domes** shield against wax and moisture but can impact sound quality if debris builds up. **In addition to daily cleaning, replace domes monthly.**

**Finding the Right Fit:** The hearing aids come with medium domes. If these don't fit, try the small or large domes included in your box to find the most comfortable size.

**Remove Dome:** Hold the hearing aid in one hand and grip the dome with the other. Twist gently back and forth while pulling to detach.



**Replace Dome:** Push the new dome securely onto the ridge at the tip of the hearing aid.



## **KEY INFO**

Always ensure the dome is fully secured before inserting the hearing aid. An improperly attached dome may detach when removing it from your ear.



**Maintenance:**

- Repairs and modifications must be performed only by authorized professionals. Any unauthorized maintenance voids the warranty.
- Do not attempt repairs yourself. For assistance or replacements, contact customer support.
- Tampering with the device or accessories is prohibited and will void the warranty.

**Storage:**

- Store devices in the charging case when not in use.
- Avoid exposing devices and the case to extreme temperatures, high humidity, or physical shocks.
- Storage environment: 14°F to 131°F, 10%-90% RH.
- Recharge devices at least every 6 months if stored for extended periods.

**Disposal:**

- Dispose of devices per the WEEE Directive 2012/19/EU. Contact your local authority for guidance on electronic waste disposal.

**Safe Usage:**

- Keep small parts out of children's reach.
- Avoid dropping, submerging in liquids, or exposing devices to strong electromagnetic fields.
- Remove devices before X-rays, MRIs, or CT scans.

- Do not use at high volume for long periods.
- Stop use if uncomfortable.
- Maintain cleanliness to prevent wax or debris buildup.
- Recycle batteries properly; do not discard them in the trash.
- Keep away from children and pets to avoid damage.
- Do not share hearing aids to prevent hearing damage or injury.
- Stop use and consult a physician if you experience any skin irritation or allergic reactions.
- Avoid strong electromagnetic fields (e.g., microwaves, ultrasonic cleaners) to prevent damage.
- Contact customer support or a hearing professional for assistance or unexpected issues.

#### **Electromagnetic Compatibility (EMC):**

- The device complies with EMC standards (EN60601-1-2) and has undergone extensive testing.

#### **Quality Assurance:**

- Manufactured under strict quality standards, the device complies with IEC60601-1, IEC60601-2-66, IEC60118-7, IEC60118-13, ISO 10993, and ISO 14971.



|                |                |
|----------------|----------------|
| Test Standard  | ANSI 3.22-2014 |
| Test Equipment | FONIX8000      |
| Test Coupler   | HA-1           |

#### **Standard test ambient conditions**

|                      |                                               |
|----------------------|-----------------------------------------------|
| Temperature          | 23°C +- 5°C (73° +- 9°F)                      |
| Relative humidity    | 0% to 80%                                     |
| Atmospheric pressure | 760 (+37.5-750) mm of Hg or 707.2 (+5-20) kPa |

#### **Technical Specifications:**

|                                 |                     |
|---------------------------------|---------------------|
| Maximum output (OSPL90)         | 107 dB              |
| Full-on gain (HFA FOG@50 dBSPL) | 29dB                |
| Total harmonic distortion       | 1%                  |
| Self-generated noise (EIN)      | 23 dB               |
| Latency                         | 1ms                 |
| Frequency response bandwidth    | f1=200Hz, f2=5600Hz |
| Frequency response smoothness   | 5dB                 |

## Symbols

|                                                                                   |                                                            |                                                                                   |                                  |
|-----------------------------------------------------------------------------------|------------------------------------------------------------|-----------------------------------------------------------------------------------|----------------------------------|
|  | Follow operating instructions                              |                                                                                   |                                  |
|  | B type applied part                                        |                                                                                   |                                  |
|  | Disposal in accordance with Directive 2012/19/EU (WEEE)    |                                                                                   |                                  |
|  | Date of manufacture                                        |  | Manufacturer information         |
|  | Serial number                                              |  | USE BY                           |
|  | Warning/Danger: Improper use might cause serious injuries  |                                                                                   |                                  |
|  | Medical Device                                             |                                                                                   |                                  |
|  | Permissible storage and transport temperature and humidity |                                                                                   |                                  |
|  | Permissible operating temperature and humidity             |                                                                                   |                                  |
|  | Keep Dry                                                   |  | Comply with Directive 2011/65/EU |
|  | Paper Packaging Recycling                                  |                                                                                   |                                  |

Notes: This hearing aid is produced under the strict quality system of GA Trading LLC.

The device fulfills provisions of IEC 60601-1, IEC60601-2-66, IEC60118-7, IEC 60118-13, ISO 10993-1/-5/-10 ISO 14971.



