

User Guide

Tianxi

OTC Hearing Aid



Please refer to this user guide before use.

FDA Warnings and Cautions



WARNING

If You are younger than 18, do not use this.

You should go to a doctor, preferably an ear-nose-throat doctor (an ENT), because your condition needs specialized care. Over-the-counter hearing aids are only for users who are 18 or older.

This OTC hearing aid is for users who are 18 and older. People who are younger than 18 with hearing loss should see a doctor, preferably an ENT, because they may need medical testing and management. Hearing loss can affect speech and learning, so professional fitting and continuing care are also important.



WARNING

When to see a doctor.

If you have any of the problems listed below, please see a doctor, preferably an ear-nose-throat doctor (an ENT).

- Your ear has a birth defect or an unusual shape. Your ear was injured or deformed in an accident.
- You saw blood, pus, or fluid coming out of your ear in the past 6 months.
- Your ear feels painful or uncomfortable.
- You have a lot of ear wax, or you think something could be in your ear.
- You get really dizzy or have a feeling of spinning or swaying (called vertigo).
- Your hearing changed suddenly in the past 6 months.
- Your hearing changes: it gets worse, then better again.
- You have worse hearing in one ear.
- You hear ringing or buzzing in only one ear.



WARNING

This hearing aid should not cause pain when inserting it.

Remove the device from your ear if it causes pain or discomfort when you insert or place it. To try again, make sure to follow the instructions. If you feel pain or discomfort again, contact the manufacturer. If your pain or discomfort does not go away, contact your hearing healthcare professional. You can also report this to FDA as an adverse event according to the instructions that appear later.



WARNING

If you remain concerned, consult a professional.

If you try this device and continue to struggle with or remain concerned about your hearing, you should consult with a hearing healthcare professional.

These precautions point out situations which could cause minor or moderate injuries to the user or to the hearing aid.

- You should familiarize yourself fully with the entire contents of the user guide before using the hearing aid to ensure personal safety and correct use.
- Contact the manufacturer for support if you experience unexpected operations or performance.
- In some acoustic environments, over amplification can cause discomfort and damage to the patient's hearing.
- A hearing aid should only be used as directed by your physician or hearing healthcare provider.
- Do not use your hearing aid in explosive or oxygen-enriched areas.
- Allergic reactions to hearing aids are unlikely. However, if you experience itching, redness, soreness, inflammation, or a burning sensation in or around your ears, inform your hearing healthcare provider and contact your physician.
- You may notice accelerated accumulation of earwax when using your hearing aid. This is normal, but see your hearing healthcare provider if you continue to have difficulties.

- In the unlikely case that any parts remain in the ear canal after the removal of the hearing aid, contact a physician immediately. Report the event to your hearing healthcare provider as soon as possible.
- Remove your aids for CT and MRI scans or other electromagnetic procedures.
- Never leave your hearing aids or batteries where small children and pets can reach them. If batteries are swallowed or inhaled, contact the local poison control center immediately.
- If you feel pain in or behind your ear, if it is inflamed, if skin irritation or accumulations of ear wax occur, please check with your hearing healthcare provider or physician.
- Your hearing aids have been adjusted for your hearing needs, please do not share your hearing device with other people.

 CAUTION

This is not hearing protection.

If you are exposed to excessively loud sounds (regardless of duration), you should remove this device. If you are in a noisy environment, you should use appropriate hearing protection earplugs instead of wearing this device. Generally speaking, if you need to use hearing protection earplugs in a noisy environment, you should remove this device and wear professional hearing protection earplugs.

 CAUTION

The sound output should not be uncomfortable or painful.

You should turn down the volume or remove the device if the sound output is uncomfortably loud or painful. If you consistently need to turn the volume down, you may need to further adjust your device.



CAUTION

You might need medical help if a piece gets stuck in your ear.

If any part of your hearing aid, like the ear plug, gets stuck in your ear, and you cannot easily remove it with your fingers, get medical help as soon as you can. You should not try to use tweezers or cotton swabs because they can push the part farther into your ear, injuring your eardrum or ear canal, possibly seriously.



Note: If you remain concerned, consult a professional.

If you try this device and continue to struggle with or remain concerned about your hearing, you should consult with a hearing healthcare professional.

❗ Note: What you might expect when you start using a hearing aid.

A hearing aid can benefit many people with hearing loss. However, you should know it will not restore normal hearing, and you may still have some difficulty hearing over noise. Further, a hearing aid will not prevent or improve a medical condition that causes hearing loss.

People who start using hearing aids sometimes need a few weeks to get used to them. Similarly, many people find that training or counseling can help them get more out of their devices.

If you have hearing loss in both ears, you might get more out of using hearing aids in both, especially in situations that make you tired from listening - for example noisy environments.

ⓘ Note: Tell FDA about injuries, malfunctions, or other adverse events.

To report a problem involving your hearing aid, you should submit information to FDA as soon as possible after the problem. FDA calls them "adverse events", and they might include: skin irritation in your ear, injury from the device (like cuts or scratches, or burns from an overheated battery), pieces of the device getting stuck in your ear, suddenly worsening hearing loss from using the device, etc.

Table of Contents

FDA Warnings and Cautions	1
Introduction	11
What's in the Box	14
Getting to Know Your Hearing Aids	15
Preparing Your Hearing Aids for Use	17
Charging Your Hearing Aids	19
Turning Your Hearing Aid On and Off	24
Placing Your Hearing Aids in Your Ears	25
Remove the Hearing Aid from Your Ears	26
Adjust Your Hearing Aid's Volume and Mode	27
Cleaning and Caring for Your Hearing Aids	29
Technical Specifications	35
Troubleshooting Guidelines	38
General Warnings and Cautions	40
Symbols	43

Introduction

Thank you for choosing this hearing aid. We recommend that you wear your hearing aid every day.

NOTE:

Please read this user guide carefully before using the hearing aid.

Intended Use

This hearing aid is designed to amplify sound for individuals aged 18 and above who perceive they have mild to moderate hearing loss.

User Profile

The hearing aid is intended to be used by adults 18 years of age or older.

Therapeutic Indications

Perceived mild to moderate hearing impairment.

Contraindications

You should see a physician if you have:


- Trauma, bleeding or visible deformity of the ear. Fluid, drainage, or pain/discomfort in one or both ears within the past 90 days.
- Recent or recurrent ear infection, or a plugged-up fullness feeling.
- Sudden or rapid progression of hearing loss within the past 90 days in one or both ears.
- Hearing loss on one side that has worsened in the past 90 days.
- Spells of acute or chronic dizziness, poor dexterity, poor vision, or significant dementia.
- Excessive ear wax build-up, or a history of excessive wax build-up, itching, or a foreign object in the ear canal.
- Sudden onset or rapid worsening of ringing in one or both ears within the past 90 days.

Side Effects


If you experience side effects, contact a hearing care professional or a physician.
Possible side effects from wearing a hearing aid may be:

- Dizziness
- Tinnitus
- Perceived worsening of hearing loss
- Nausea
- Skin reaction
- Ear wax accumulation

What's in the Box

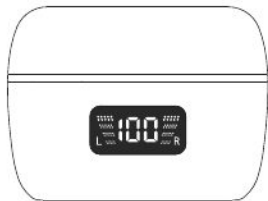

Charging Case

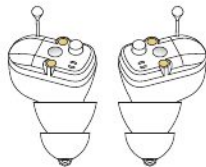


Hearing Aid



Domes




Cleaning Brush

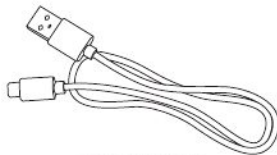


Charging Line

Getting to Know Your Hearing Aids

Hearing Aids

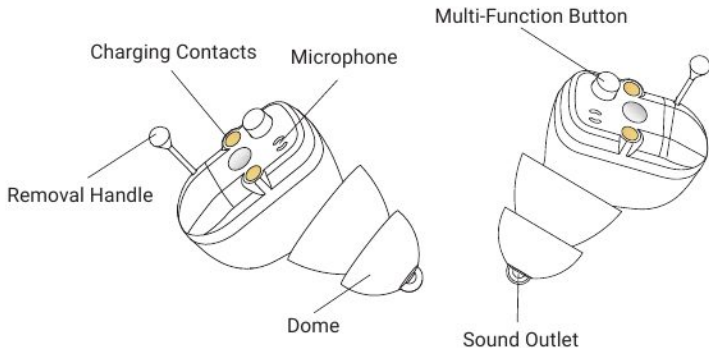

Charging Case

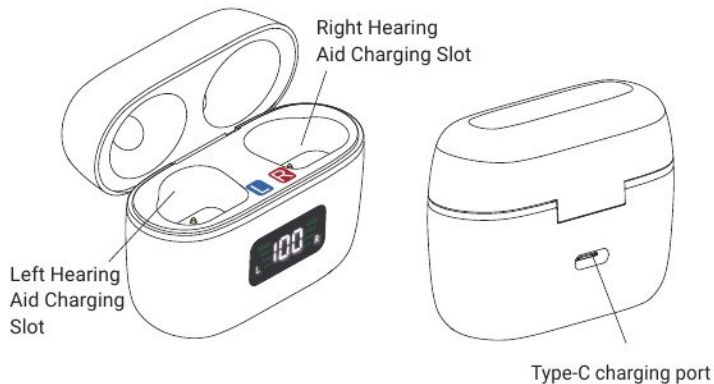

Preparing Your Hearing Aids for Use



Battery Warnings

- This hearing aid is powered by a non-removable rechargeable lithium-ion battery.
- Contact with or ingestion of the battery may cause serious injury.
- Never put the rechargeable hearing aid in your mouth.
- If the housing of the rechargeable hearing aid is damaged, the internal rechargeable battery may leak. In such cases, do not attempt to use the hearing aid.
- Battery leakage may cause chemical burns. If you come into contact with battery leakage, rinse immediately with warm water. If chemical burns, redness, or skin irritation occur due to battery leakage, seek medical assistance.
- Batteries are harmful to the environment. Therefore, never attempt to incinerate the battery. Dispose of used rechargeable hearing aids in accordance with the regulations of your country or region.

- For safety reasons, please use the charging cable supplied with the hearing aid and a standard-compliant 5V/1A charger.
- When traveling, please consult the airline about the transportation regulations for hearing aids and chargers, as the charger of your hearing aid contains a lithium-ion battery.

NOTE:


When not using the hearing aids, you can place them back in the charging case to turn them off and charge them.

Charging Your Hearing Aids

We recommend that you fully charge the hearing aid before using it. This is merely a precautionary measure to ensure you do not accidentally run out of power.


NOTE:

1. Hearing aids leave the factory with an insulating sticker on their charging contacts to prevent battery pre-discharge and preserve factory charge. Remove the sticker before first use.

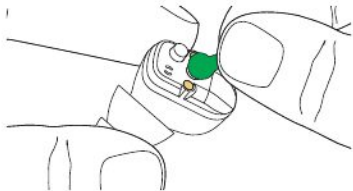

2. Before putting your hearing aids into the charger, ensure that they are dry.
3. When charging, place the hearing aid and its charging case in a dry, well-ventilated and cool area. Do not charge in areas with direct sunlight, high temperature, humidity or heavy dust to prevent the device from getting damp, parts aging or short circuits.

Charging Indicator Light

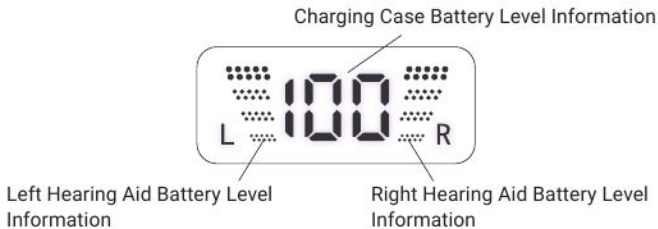

Number of lights				
Battery level	0 – 25%	26 – 50%	51 – 75%	76 – 100%

Charging the Hearing Aid Using the Charging Case

The charging case has a power storage function. When the charging case has power, simply place the hearing aid into the case correctly, and the hearing aid will start charging automatically. During charging, the hearing aid's battery indicator light flashes green; once fully charged, the indicator light turns off.




Connecting Power to Charge the Charging Case and Hearing Aid

When the charging case runs out of power, you can use an external power source to charge the charging case. When the charging case is connected to an external power source, the hearing aids can also be placed in the charging case to charge together.

Step 1: Connect the charging cable to the power adapter and the charging case.

Step 2: Plug the adapter into a power outlet.
This starts charging the case (and the hearing aid, if placed inside).

Step 3: During the case's charging, its power indicator number rises to 100—this means charging is complete.

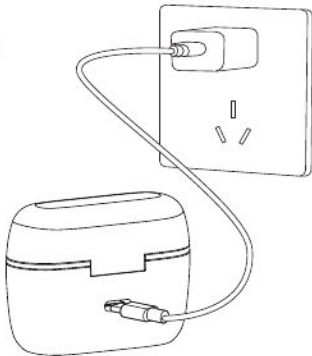

Charging Information

Hearing Aid Battery Capacity	18mAh
Hearing Aid Charging Time	Approximately 2 hour
Hearing Aid Battery Usage Time	>18 hours
Charging Case Battery Capacity	320mAh
Charging Case Charging Time	Approximately 2 hour
Charging Times of a Fully Charged Case	Can charge a pair of hearing aids 5 times

Turning Your Hearing Aid On and Off

Turning On

Remove the hearing aids from the charging case; they power on automatically after a prompt tone.

Turning Off

Place your hearing aid back into the charging case, and it will turn off automatically and start charging.



Placing Your Hearing Aids in Your Ears

Step 1: Distinguish between the left and right hearing aids by checking the "L/R" markings on the device body (L stands for the left ear, R stands for the right ear).

Step 2: Hold the hearing aid and insert it into the ear canal.

Step 3: Fine-tune the position of the hearing aid until it fits comfortably and securely in the ear without causing pressure.

NOTE:

When worn correctly, the multi-function button is located above the microphone.



Remove the Hearing Aid from Your Ears

Gently pull the removal handle to take the hearing aid out of your ear.



Adjust Your Hearing Aid's Volume and Mode


Adjust Volume

Short-press the multi-function button to adjust volume cyclically (9 levels in total; Level 1 is minimum, Level 9 maximum):

Each short press raises the volume by one level accompanied with a single "beep" prompt tone. Upon a double "beep-beep" alert, the volume reaches maximum.

Another short press at maximum volume cycles the volume back to the minimum, and the adjustment repeats in loop.

The volume follows a looping sequence: 1-2-3...9-1-2-3...9-1-2...

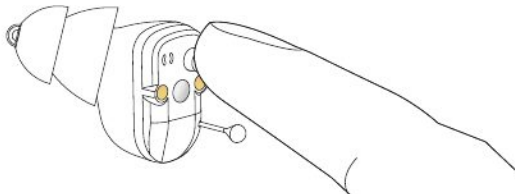

Adjust Mode

The hearing aid comes with 4 preset modes (Mode 1: Quiet Mode, Mode 2: Noise Reduction Mode, Mode 3: Outdoor Mode, Mode 4: Music Mode). Press and hold the "Multi-function Button" for about 2 seconds; after hearing a prompt tone, the hearing aid will switch modes in a loop. Different prompt tones correspond to different modes (as shown in the table below):

Prompt Beep	Corresponding Mode	Applicable Environment
1 beep	Mode 1 (Quiet Mode)	Home, 1-on-1 office talks ...
2 beeps	Mode 2 (Noise Reduction Mode)	Restaurants, cafes ...
3 beeps	Mode 3 (Outdoor Mode)	Busy streets, outdoor markets ...
4 beeps	Mode 4 (Music Mode)	Concert halls, music venues...

Cleaning and Caring for Your Hearing Aids

General Instructions for Care and Maintenance

To ensure you get the highest quality experience and longest useful lifetime out of your hearing aids, it is important to clean and care for them.

NOTE:

- After removing the hearing aid, wipe it with a soft cloth to keep it clean and dry.
- Apply cosmetics, perfume, aftershave, hairspray, moisturizer, and other such products first, then put on the hearing aid. These products may damage the hearing aid or cause discoloration.
- Do not immerse the hearing aid in liquid.
- Keep the hearing aid away from excessive heat and direct sunlight.
- Avoid contact with liquids; do not swim, shower, or use a sauna while wearing the hearing aid.




Warning:


- Do not immerse the hearing aid in liquid.
- Keep the hearing aid away from excessive heat and direct sunlight.
- Avoid contact with liquids; do not swim, shower, or use a sauna while wearing the hearing aid.

Daily Care and Maintenance

It is important to keep your hearing aids clean and dry on a daily basis. Use a cloth and brush to clean the hearing aids.

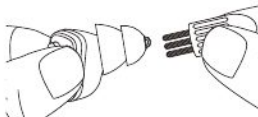

Clean the dome

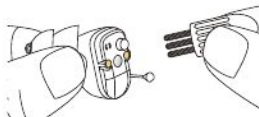


Clean the microphone port and
charging contacts

Changing Your Hearing Aid Domes

You should replace the dome every 3 months or more frequently. This can effectively prevent earwax and sweat from adhering to the dome and breeding bacteria, reduce the risk of ear canal infections, prevent the dome from aging and deforming, ensure a proper fit to avoid sound leakage and whistling, and guarantee clear hearing assistance.

Step 1: Grip the body of the hearing aid and remove the dome.

Step 2: Select ear plug that fit your ear canal and insert them into their original positions.

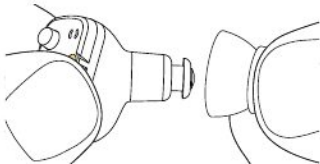



NOTE:

- Incorrect replacement of the dome may cause the dome to be left in the ear when removing the hearing aid.
- If you suspect that you have a dome, wax filter, or other object in your ear canal, consult a hearing care professional or a physician. These objects can be harmful and can cause an infection in your ear.

Technical Specifications

MAX Output (OSPL90)	<117 dB
HF Average Output (OSPL90)	102±4 dB
Full-on Gain (FOG50)	NA
HFA Full-on Gain (FOG50)	29±5 dB
Equivalent Input Noise	≤32 dB
Frequency Range	F1<250 Hz F2>5000 Hz
Total Harmonic Distortion	≤3%

Classification:

The hearing aids are classified as a Type B Applied Part, IEC 60601-1 .



Transportation and Storage Conditions



Temperature

0 to +38 °C



Relative humidity

0 to 80 %

Operation Conditions



Temperature

-10 to +40 °C



Relative humidity

0 to 80 %



Atmospheric pressure

860 to 1060 hPa

Expected service Lifetime

The hearing aids are designed with an expected service life of 5 years (3 years for the lithium battery). They are safe for use within this period.



Waste Disposal

Charging cases, hearing aids, and charging cables should not be discarded as household waste and should be sent to appropriate collection facilities for recycling. Everything else is normal waste. Proper handling and recycling can help protect natural resources, human health, and the environment.

Troubleshooting Guidelines

Malfunction	Possible Reasons for	The Solution
No sound	Low volume	Increase the volume
	Equipment is out of power	Charge your hearing aid
	Microphone clogging	Clean microphone
	Wax guard is blocked	Replace wax guard
The voice was low, weak, or unclear	The battery is low	Charge your hearing aid
	The dome are too loose	Choose the right size to wear
	Dome or wax guard have wax on them	Clean dirt from ears and dome, Replace New wax guard
	Hearing change	Contact your doctor

Malfunction	Possible Reasons for	The Solution
Hearing aid whistles in the charging case	The hearing aid is placed out of position	Remove the hearing aid from the charging case and place the hearing aid in the charging case according to the "L,R" mark
	Low battery of the charging case	Charge the charging case
Hearing aid is in the charging case, the hearing aid charging indicator light is not on	Poor contact of charging contact	Insert the hearing aid back into the charging base. If the issue persists, ensure that the charging contacts are clean. Use the provided brush from the accessories kit to thoroughly remove any debris
	Green insulation sticker cover has not been removed	Remove and discard green insulation sticker cover

General Warnings and Cautions



Warning:

- No modification of this hearing aid is allowed.
- Seek immediate medical help if a hearing aid, any of its parts, a battery, or a magnet is swallowed, as it can cause choking and harm your health.
- Keep hearing aids, their parts, batteries, and magnets away from pets, children, and people with cognitive, intellectual, or mental health challenges.
- Always supervise children or individuals with cognitive, intellectual, or mental health challenges when they are using their hearing aids. Hearing aids contain small pieces that can be dangerous if swallowed.
- Do not wear your hearing aids while being exposed to radiation. Some types of radiation, such as from MRI or CT scanners, can affect the settings in your hearing aids, causing malfunction and potentially damage to your hearing.

- Other types of radiation, such as burglar alarms, room surveillance systems, mobile phones, metal detectors, and radio equipment will not damage your hearing aids. However, they may briefly affect the sound quality in your hearing aids and may create undesired sounds.
- Never use your hearing aids in places with explosive gases such as mines, oil fields, or similar unless these areas are certified for hearing aid use. Using your hearing aids in places that are not certified for hearing aid use can be dangerous.
- Do not dry your hearing aids in an oven, microwave oven, or other heating equipment. This will cause them to melt and may cause burns to your skin.
- In general, exposure to loud sounds can damage your hearing. This could be loud music or loud environments. You can best protect your hearing by reducing exposure to loud sounds and environments or by using hearing protection.

**Caution:**

- Use as professionally advised; improper operation may easily damage hearing.
- Do not use damaged/modified devices to prevent malfunction, hearing damage and skin scratches.
- Do not alter the shape of the device or its accessories to avoid skin reactions or scratches.
- If there is a wound on the ear/head area in contact with the device, stop using it immediately and consult a professional promptly.
- For personal use only; do not lend to others to prevent harm to their hearing.

Symbols



WARNING: Points out a situation that could lead to serious injuries.



Legal manufacturer.



Unique Device Identification.



Follow instructions for use.



Date of manufacture.



According to IEC 60601-1 the applied parts of this device are classified as type B.



Xiamen Vienatone Medical Technology Co., Ltd.

8th Floor, Block A, Jianye Building, No. 96 Xiangxing
Road, Xiamen Torch High-Tech Zone (Xiang'an)
Industrial Zone, 361101 Xiamen, Fujian, P.R. China

www.vienatone.com